


# Optical Tune-Ups

Advertisement

FOR ALL INTRAORAL CAMERAS  
REPAIRS SPOTTY AND FOGGY IMAGES



Restores your IntraOral Camera back to it's original great image

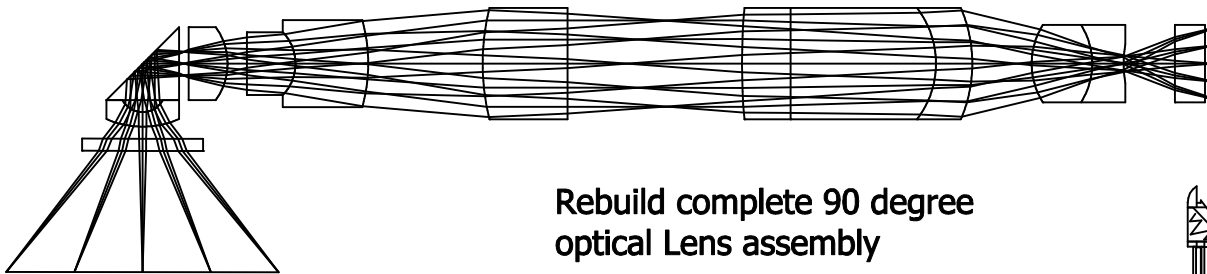


90 degree  
IntraOral Lens

## Rebuild complete 90 degree optical Lens assembly

Disassemble Optical components  
Micro clean all components  
Reassemble in a Clean Room envirement  
Recure cure bonded Optics in a UV oven  
Purge Optical assembly with Gas  
Set Focus to Manufactures Specifications  
All Lens are put through Q.C. and tested  
before shipping

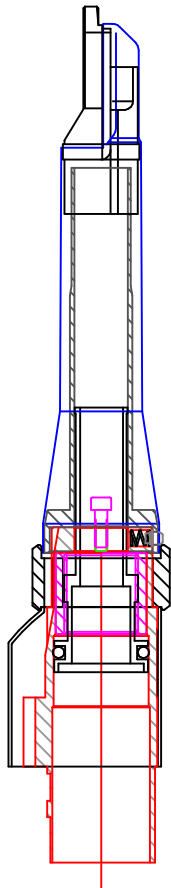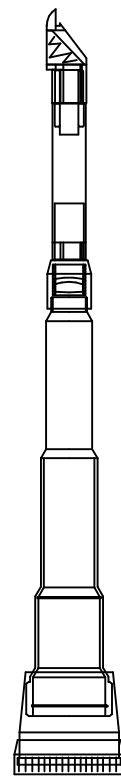

## Replace Filter Lens located in Handpiece

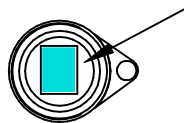

Remove faulty Filter Lens  
Micro Clean Assembly  
Purge CCD area with Gas  
Install new Filter Lens  
Inspect with Microscope for debre  
Test system for Clarity

ALL OPTICAL REPAIRS ARE PERFORMED  
IN A CLEAN ROOM ENVIRONMENT

Manufactured by

**IMAGIN SYTEMS CORP.**

815 CHERRY LANE, SAN CARLOS, CALIFORNIA 94070, TEL. (650) 596-1001

9-1-04